

LOMBA LS FIX

Back support

SRT 137

Indications

Used for incipient intervertebral disc herniation, sciatica, moderate low back pain, post-operative, osteoporosis arthritis and moderate stabilisation of the lumbar spine.

Description

This lumbosacral brace is made with breathable elastic fabric, provided with:

- 4 removable, malleable contoured metal stays located on the back, to support the lumbar area
- plastic reinforcement on the front section provides counter pressure to support the lumbar area
- 25 cm elastic straps around the waist provide adjustable compression and stability


Hook and loop fastener on the front side, with finger loops for easy application.

Width: front 19 cm, back 30 cm.

4 metal stays stabilize the lumbar region, slightly reducing lordosis, and the compression exerted by the elastic bands reduce muscle tension and pain.

Composition


Breathable elastic fabric: 50% polyester, 33% polyamide, 17% elastane, designed to allow high air permeability in both directions.

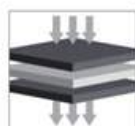

Air Transfer



- Available in Black

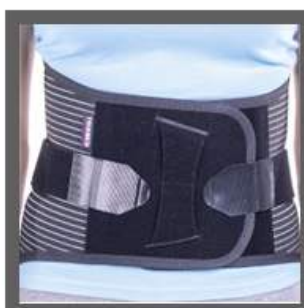


- Finger loops for easy application

- Elastic Straps are adjustable for enhanced support and compression



- Rigid reinforcement to the front for increased support

Sizes	1	2	3	4	5	6
Waist circumference (cm)	65-75	75-85	85-95	95-105	105-115	115-125
Reference	SRT137M1	SRT137M2	SRT137M3	SRT137M4	SRT137M5	SRT137M6

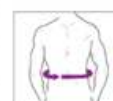

*Measurement taken from the circumference of the waist