

STOMATEX

Abdominal stoma support

SRT 106

Indications

Used after abdominal surgery for people with colostomy or urostomy. Used to prevent and treat parastomal hernia.

Description

The support is made of elastic fabric and firm textile. The frontal, firm textile part is fitted with an opening for the colostomy bag, made with a 25 mm wide and 75 mm inside diameter polyethylene ring.

Hook and loop fastener on the front side. Width: 27,5 cm. The firm anterior panel supports the abdominal muscle around the stoma and can be used to prevent or treat a parastomal hernia. It is a semirigid traditional corset that due to the elastic component, allows normal abdominal breathing.


Composition

Elastic fabric: 50% polyester, 33% polyamide, 17% elastane, with a structure that allows high air permeability in both directions.



Air Transfer

The triple-layer laminated textile construction allows for moisture regulation, ensuring a quick transfer of moisture away from the body.

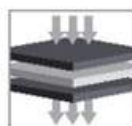

Air Transfer



Moisture Transfer



- Available in Black or Beige




- traditional semi-rigid corset ensures uniformed compression on the entire covered surface
- the bag can be easily attached and accessed



- supportive firm textile
- elastic fabric for a snug fit while not impeding function

Sizes	1	2	3	4	5	6
Waist circumference (cm)	65-75	75-85	85-95	95-105	105-115	115-125
Reference	SRT106M1	SRT106M2	SRT106M3	SRT106M4	SRT106M5	SRT106M6

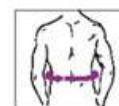

*Measurement taken from the circumference of the waist.