

TRIASOR 02

Back support **SRT 108**

Indications

Used for incipient intervertebral disc herniation, moderate low back pain, postoperative, and for osteoporosis, sciatica and arthritis.

Description

The support is made of breathable textile, and incorporates:

- 4 integrated plastic stays for support
- one plastic lumbar pad with cotton liner
- 2 6 cm criss-cross elastic straps around the waist providing increased stability and compression.

Hook and loop fastener on the front side. Width: front 19 cm, back 30 cm. The inside pad has 8 leaf springs arranged on 2 rows that have a lumbar massage effect.

4 plastic stays stabilize the lumbar region, slightly reducing lordosis, and the compression exerted by the elastic bands reduces muscle tension and pain.

Composition

Breathable three-layer laminated moisture-wicking textile.

50% polyester, 33% polyamide, 17% elastane






Air Transfer Water Transfer



Available in Black



Reference	SRT108M1	SRT108M2	SRT108M3	SRT108M4	SRT108M5	SRT108M6
Waist circumference (cm)	65-75	75-85	85-95	95-105	105-115	115-125
Sizes	1	2	3	4	5	6

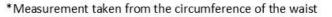





- elastic straps around the waist providing a better stability and compression
- the inside pad has leaf springs that produce a lumbar massage effect



 Supportive lumbar pad is fitted with springs which produce a lumbar massage effect.

